


Mammography X-Ray QA/QC

# **Technical Accessories**

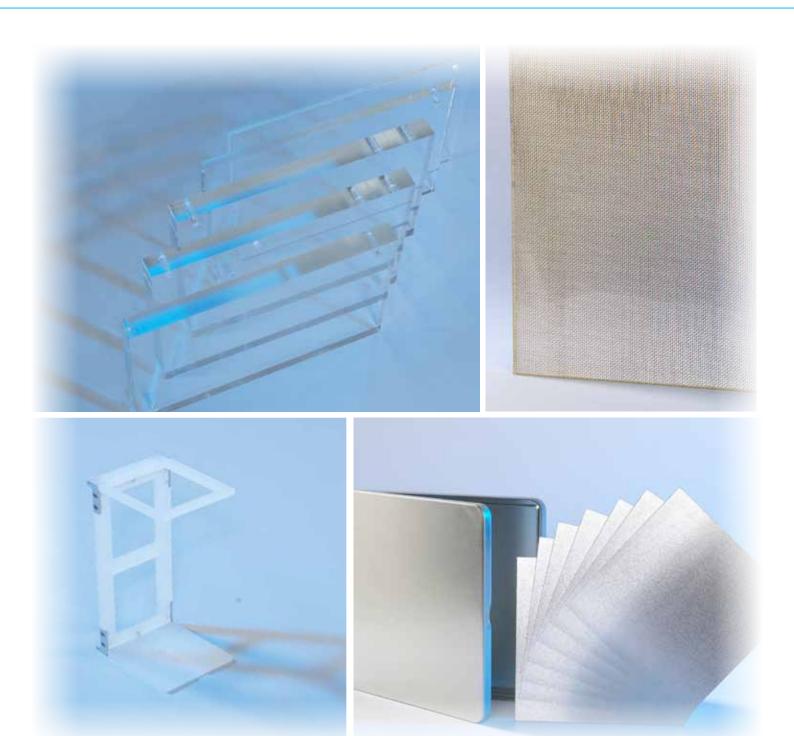

**深圳市展业达鸿科技有限公司** 吕先生:15920060912(微信同号) 0755-22934005(座机) 地址:深圳市福田区八卦二路八卦岭工业区615栋419 邮箱:hongqi@thingstet.com 网址:www.thingstest.com Technical Accessories

AEC Test Set / HVL Stand / HVL Set

## AEC Test Set

The Mammography AEC Test Set consists of 5 PMMA plates of varying thickness. The set is used to test the automatic exposure control features of mammography equipment according IEC 61223-3-1.

For checks on the AEC feature, test setups with attenuation ranging from 5 - 75 mm can be realised. Additional slabs to add to the set are available on request. The outer dimensions of the plates are  $180 \times 240$  mm.

**Delivery includes** 

- ✓ 3 Slabs 180 x 240 x 20 mm
- ✓ 1 Slab 180 x 240 x 10 mm
- ✓ 1 Slab 180 x 240 x 5 mm
- ✓ Manual
- ✓ Transport Case

### **HVL** Stand

The HVL Stand is designed to make HVL measurements straightforward. The collapsible stand is easy to transport and fast to set up. It allows a reproducible and secure positioning of filters and dosemeter detector.

**Delivery includes** 

- ✓ HVL Stand
- ✓ Manual
- ✓ Soft Bag

## HVL Set

The HVL Set consists of 8 filters with varying thicknesses: 1 x 2.0 mm / 2 x 1.0 mm / 2 x 0.5 mm / 1 x 0.3 mm / 2 x 0.1 mm = overall 5.5 mm. The filters have a guaranteed purity of 99.5% which ensures accurate measurements. The format of all filters of the set is 10 x 10cm. They are fully compatible with the filter insert of the QUART HVL Stand.


**Delivery includes** 

- ✓ HVL Filter Set
- ✓ Manual
- ✓ Filter Case







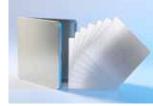



#### 深圳市展业达鸿科技有限公司

吕先生:15920060912(微信同号) 0755-22934005(座机) 地址:深圳市福田区八卦二路八卦岭工业区615栋419 邮箱:hongqi@thingstet.com 网址:www.thingstest.com

© 2007 QUART GmbH

## **Technical Accessories**

Pinhole Camera / Focal Spot Test Tool / Cassette Test Tool

## Pinhole Camera (mammotest 1.0)

The QUART mammotest 1.0 pinhole camera is designed to verify shape and size of focal spots of a mammography x-ray unit >0.1 mm. To acquire the image data, no modification or mechanical intervention in the equiment is necessary!

For checks on the focal spot properties, the mammotest 1.0 is more suitable than a star pattern or slit camera. Our tool is capable of imaging the actual shape of the focal spot in one exposure, including hot spots and image areas with less intensity.

**Delivery includes** 

✓ QUART mammotest 1.0 Pinhole Camera

✓ Manual

## Focal Spot Test Tool (mammotest 2.0)

The QUART mammotest 2.0 focal spot test tool is designed to visualise the focal spot distribution within a digital or conventional mammography image.

The check mainly verifies, if and to what extent does the position of the focal spot affect the equipment's imaging properties in the corresponding region of interest.

**Delivery includes** 

- ✓ QUART mammotest 2.0 Focal Spot Test Tool
- ✓ Manual

## **Cassette Test Tool**

The QUART KAT 24x30 is used to check the film contact at conventional mammography equipment with film/screen technology. Flaws would be revealed by variations in density on the developed film, i.e. showing light or dark areas.

**Delivery includes** 


✓ QUART KAT 24x30 Film Contact Test Tool

✓ Manual

Also available QUART mam/digi QUART dido/time M QUART MKmam 18x24 QUART MKmam 24x30 QUART al20

Phantom for QA/QC in Mammography Constancy Test Dosemeter Measuring Cassette for Dosemeter Detectors Measuring Cassette for Dosemeter Detectors Standard Compliant Added Filtration



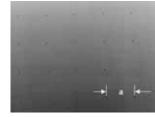





© 2007 QUART GmbH



深圳市展业达鸿科技有限公司

吕先生:15920060912(微信同号) 0755-22934005(座机) 地址:深圳市福田区八卦二路八卦岭工业区615栋419 邮箱:hongqi@thingstet.com 网址:www.thingstest.com









QUART GmbH • Kirchenweg 7 • D-85604 Phone +49 (0) 8106 / 249 118 • Fax +49 (0) 81 深圳市展业达鸿科技有限公司

吕先生:15920060912(微信同号) 0755-22934005(座机) 地址:深圳市福田区八卦二路八卦岭工业区615栋419 邮箱:Nongqi@thingstet.com 网址:www.thingstest.com