

# D100 Water Tank

**QRM-Water-Tank provides the opportunity to evaluate noise and homogeneity in CT images. The water tank was designed to fit in the standard thorax and abdomen phantoms. It can be used as well to measure sample probes in water.**

The Phantom can be attached to the standard QRM-Thorax and QRM-Abdomen, but can as well be used as a stand alone phantom.


The QRM-Water-Tank can be filled with water and offers direct access to place smaller sample probes in the cylindrical unit.

It provides the ability to evaluate noise and homogeneity of a 3D scan.


The phantom is as well part of our motion simulators: Sim4D and Sim2D!

## Specifications


Cylinder diameter .....	100 mm
Phantom length .....	228 mm
Phantom weight .....	0.6 kg
Cylinder wall thickness .....	3.0 mm
Material .....	PMMA

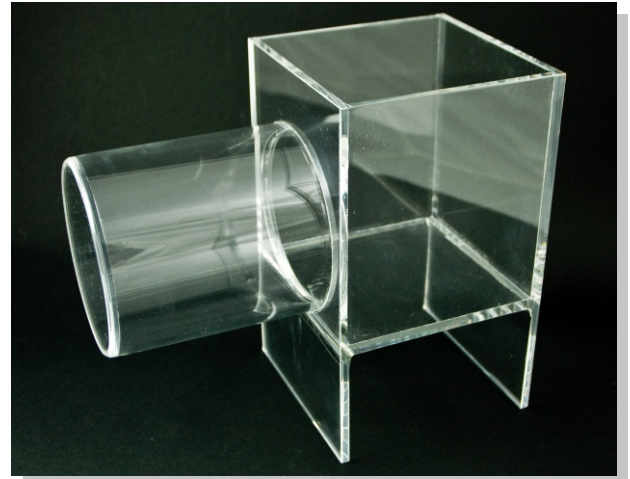

*QRM-Water-Tank*

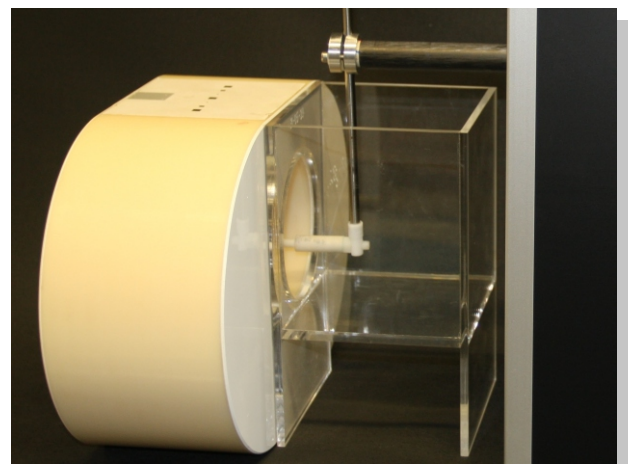

*Water-Tank attached to QRM-Thorax*

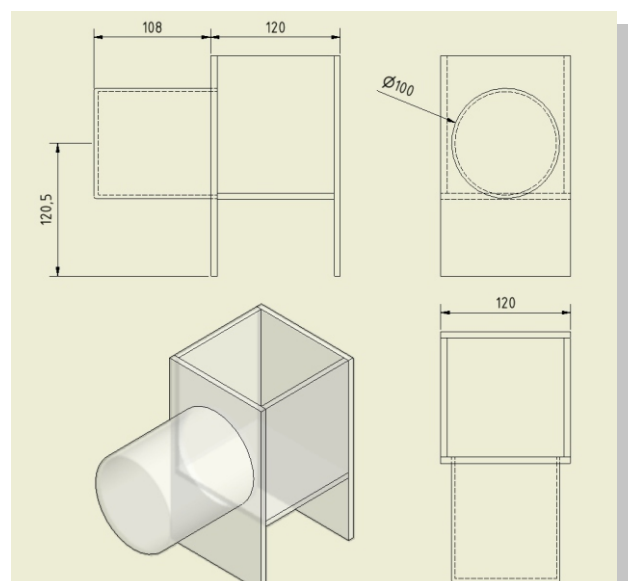

**深圳市展业达鸿科技有限公司**

Line sketch of the Water Tank

地址: 深圳市福田区八卦二路八卦岭工业区615栋419

邮箱: hongqi@thingstet.com

网址: www.thingstet.com